



How Gentle hair removal treatments work

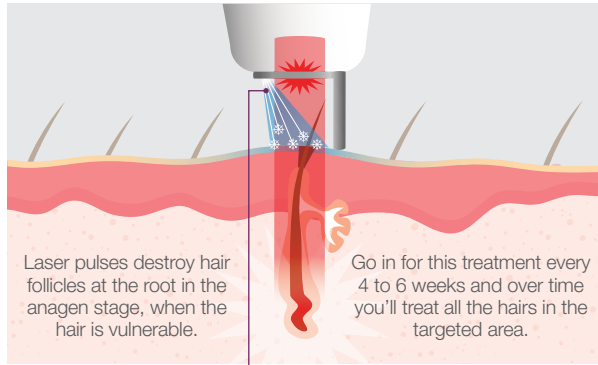

The Gentle system uses a patented Dynamic Cooling Device™ (DCD™), ensuring treatments are as comfortable as they are effective. This sprays the upper layers of the skin with a cooling burst of cryogen milliseconds before and after the laser pulse⁴.

The Gentle system is a world leader, with capabilities treating beyond laser hair removal

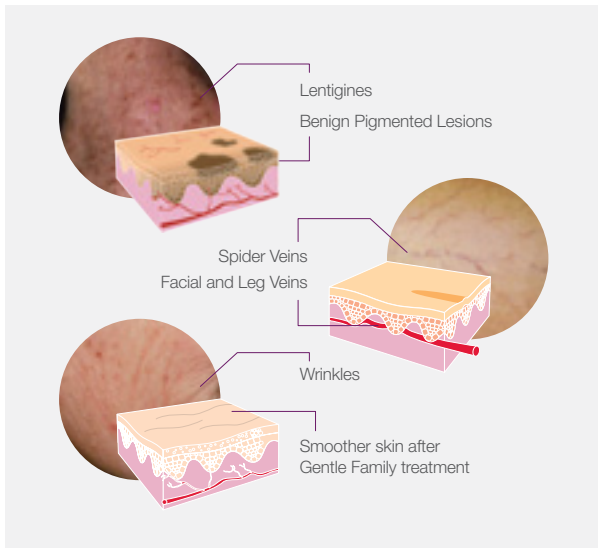

Make sure your laser hair removal treatments are Gentle

**Schedule a personal consultation
with your practitioner today!**

For more information about Gentle visit
candelamedical.com

References: 1. GentleMax Pro Plus FDA 510(k) K201111, May 2014. 2. Lee, S.J., 2021. An optimal treatment for facial pigmented lesions, Candela clinical insight. 3. Lee, J. H., et al. 2015. Efficacy of the long-pulsed 1064-nm neodymium:yttrium-aluminum-garnet laser (LPND) (rejuvenation mode) in the treatment of papulopustular rosacea (PPR): A pilot study of clinical outcomes and patient satisfaction in 30 cases, J. AM. ACAD. Dermatology, 73(2), pp. 333-336. 4. DCD cooling allows for higher fluences to be used for permanent hair reduction (1999 Dynamic Cooling Device Protects Epidermis for Effective Laser Treatment, Stephen Eubanks, MD Dermatology & Laser Center, Denver, CO. 5. Proprietary Syneron Candela survey conducted by Wizer via Nielsen's online physicians' panel, among dermatologists, plastic surgeons and medspas, 2015. 6. Eremia S, Li C, Newman N. Dermatol Surg. 2001 Nov Laser hair removal with alexandrite versus diode laser using four treatment sessions: 1-year results 7. Kim YK, Kim DY, Lee SJ, et al. J Eur Acad Dermatol Venerol. 2014 Aug. Therapeutic efficacy of long-pulsed 755-nm alexandrite laser for seborrheic keratoses 8. Lee JH, Kim M, Bae JM, et al. J Am Acad Dermatol. 2015 Aug Efficacy of the long-pulsed 1064-nm neodymium:yttrium-aluminum-garnet laser (LPND) (rejuvenation mode) in the treatment of papulopustular rosacea (PPR): A pilot study of clinical outcomes and patient satisfaction in 30 cases 9. Lee JH1, Kim M1, Bae JM1, Cho BK1, Park HJ2., Efficacy of the long-pulsed 1064-nm neodymium:yttrium-aluminum-garnet laser (LPND) (rejuvenation mode) in the treatment of papulopustular rosacea (PPR): A pilot study of clinical outcomes and patient satisfaction in 30 cases. J Am Acad Dermatol. 2015 Aug;73(2):333-6. doi: 10.1016/j.jaad.2015.05.030 10. Paranjape A, Schomacker K. Executive Summary for 2ms vs 3ms Pulse Width for Gmax PRO, December 2018

Disclaimer: All contents of this material are for informational purposes only and provided by Candela without warranties of any kind. Healthcare professionals are solely responsible for making their own independent evaluation as to the suitability of any product for any particular purpose and in accordance with country specific regulations. The availability of products and the indications mentioned in this material is subject to the regulatory requirements and product registration status in each country. Refer to the User Manual for country specific indications. Products and technical specifications may change without notice. Please contact Candela for more details.

© 2022 Candela Corporation. This material contains registered and unregistered trademarks, trade-names, service marks and brand names of Candela Corporation and its affiliates. All other trademarks are the property of their respective owners. All rights reserved.



GentleMax Pro Plus™
Excellence, evolved.

**Fast, proven, permanent
hair reduction and so much more**

GentleMax Pro Plus™

The GentleMax Pro Plus™ system is the latest evolution of the Gentle family lasers with increased capabilities on an evolved platform. It consists of the same dual-wavelengths of both 755 nm Alexandrite and the 1064 nm Nd:YAG as its predecessor and has further expanded in its ability to treat more efficiently, comfortably and can even successfully treat fine and residual hair with clinically proven results^{1-3, 7-10}.

What are the treatments like?

The laser delivers an intense beam of light, which is absorbed by the skin and converted to heat energy. Some patients may feel a slight discomfort with the laser pulse, also described as a snap of a rubber band on the skin.

How many treatments will I require?

The number of sessions depends on the desired treatment area and the condition of your skin.



Wrinkle Reduction



Treatment of Benign Pigmented Lesions



Permanent Hair Reduction
(Including fine hair¹⁰)



Treatment of Vascular Lesions



Treatment of Onychomycosis
(Nail Fungus)

What precautions should I take after treatment?

After your laser treatment, prevent sun exposure by applying sunblock (SPF 30 or above). In the first few days after treatment, don't scrub or rinse the treated areas with abrasive skin cleansers. Your practitioner may advise additional precautions.

Can all skin types be treated?

Gentle offers effective laser hair removal treatments on the widest variety of skin and hair types. Your practitioner will advise on the number of treatments for permanent hair reduction, while improving skin quality and revealing silky, smooth skin.



Laser Hair Removal



Fine Hair Removal



Benign Pigmented Lesions



Pigmented Lesion & Vascular Lesion Removal

